

Evaluation of a COVID-19 IgM and IgG rapid test; an efficient tool for assessment of past exposure to SARS-CoV-2

Tove Hoffman, Karolina Nissen, Janina Krambrich, Bengt Rönnberg, Dario Akaberi, Mouna Esmaeilzadeh, Erik Salaneck, Johanna Lindahl & Åke Lundkvist

To cite this article: Tove Hoffman, Karolina Nissen, Janina Krambrich, Bengt Rönnberg, Dario Akaberi, Mouna Esmaeilzadeh, Erik Salaneck, Johanna Lindahl & Åke Lundkvist (2020) Evaluation of a COVID-19 IgM and IgG rapid test; an efficient tool for assessment of past exposure to SARS-CoV-2, *Infection Ecology & Epidemiology*, 10:1, 1754538, DOI: [10.1080/20008686.2020.1754538](https://doi.org/10.1080/20008686.2020.1754538)

To link to this article: <https://doi.org/10.1080/20008686.2020.1754538>



© 2020 The Author(s). Published by Informa UK Limited, trading as Taylor & Francis Group.



Published online: 14 Apr 2020.



[Submit your article to this journal](#)



Article views: 9318



[View related articles](#)









[View Crossmark data](#)

CLINICAL RESEARCH ARTICLE



Evaluation of a COVID-19 IgM and IgG rapid test; an efficient tool for assessment of past exposure to SARS-CoV-2

Tove Hoffman ^a, Karolina Nissen^b, Janina Krambrich ^a, Bengt Rönnberg ^{a,c}, Dario Akaberi^a, Mouna Esmaeilzadeh^a, Erik Salaneck ^{a,b}, Johanna Lindahl ^{a,d,e} and Åke Lundkvist ^a

^aDepartment of Medical Biochemistry and Microbiology, Zoonosis Science Center (ZSC), Uppsala University, Uppsala, Sweden;

^bDepartment of Medical Sciences, Infectious Diseases Uppsala University, Uppsala, Sweden; ^cLaboratory of Clinical Microbiology,

Uppsala University Hospital, Uppsala, Sweden; ^dDepartment of Clinical Sciences, Swedish University of Agricultural Sciences, Uppsala,

Sweden; ^eDepartment of Biosciences, International Livestock Research Institute, Hanoi, Vietnam

ABSTRACT

COVID-19 is the most rapidly growing pandemic in modern time, and the need for serological testing is most urgent. Although the diagnostics of acute patients by RT-PCR is both efficient and specific, we are also crucially in need of serological tools for investigating antibody responses and assessing individual and potential herd immunity. We evaluated a commercially available test developed for rapid (within 15 minutes) detection of SARS-CoV-2-specific IgM and IgG by 29 PCR-confirmed COVID-19 cases and 124 negative controls. The results revealed a sensitivity of 69% and 93.1% for IgM and IgG, respectively, based solely on PCR-positivity due to the absence of a serological gold standard. The assay specificities were shown to be 100% for IgM and 99.2% for IgG. This indicates that the test is suitable for assessing previous virus exposure, although negative results may be unreliable during the first weeks after infection. More detailed studies on antibody responses during and post infection are urgently needed.

ARTICLE HISTORY

Received 6 April 2020

Accepted 7 April 2020

KEYWORDS

COVID-19; sars-CoV-2; rapid test; IgM; IgG; diagnostics

Background

In late 2019, a cluster of cases of viral pneumonia of unknown aetiology was reported in Wuhan, Hubei Province, China. This new viral pneumonia, COVID-19 (Coronavirus Disease 2019), caused by the novel SARS-CoV-2 (Severe Acute Respiratory Syndrome Coronavirus-2), spread rapidly and developed into a global pandemic within three months from its initial detection [1–3]. Among other symptoms, those of COVID-19 often include fever and dry cough, which resemble respiratory illnesses caused by other viruses or bacteria [4–7]. Due to the overlapping manifestations, clinical diagnosis becomes problematic, especially during seasonal flu [8], why confirmation of COVID-19 depends on the detection of SARS-CoV-2 nucleic acid by reverse-transcriptase polymerase chain reaction (RT-PCR).

More than 1.26 million cases of COVID-19 in > 200 countries and territories, with more than 66.000 human deaths, have been reported ([9], 5 April 2020). Due to the limited testing in many geographical regions, it is clear that the total number of actual COVID-19 cases is much higher than the number of confirmed ones. In most of the confirmed COVID-19 cases, the patients are symptomatic showing fever, dry cough, and pneumonia, but also more

atypical symptoms such as gastrointestinal manifestations as well as anosmia and ageusia. However, the virus has been detected in completely asymptomatic individuals, e.g. in a recent study from Italy, showing that 44% of the laboratory-confirmed cases lacked symptoms [10]. The knowledge concerning the actual number of asymptomatic vs. symptomatic infections is still limited. The same is true for the potentially growing herd immunity, where almost no data is available to date.

In the present study, we evaluated a commercially available assay, the COVID-19 IgG/IgM Rapid Test Cassette (Zhejiang Orient Gene Biotech Co Ltd, Huzhou, Zhejiang, China), developed for detection of SARS-CoV-2-specific antibodies.

Material and methods

Serum samples

Capillary blood samples or serum from 29 PCR-confirmed COVID-19 patients or convalescents, and capillary blood samples from 24 healthy volunteers, without any known history of SARS-CoV-2 infection/COVID-19, were included in the study. Anonymous blood donor sera from healthy adults ($n = 80$) and 20 serum samples from babies (6–12 months) collected

before or during 2018 from the Uppsala Academic Hospital were used as negative controls. Clinical samples that had been deposited in Uppsala Biobank were anonymized and used in accordance with local ethical guidelines. They were all used with informed consent according to the Swedish Biobank law, which allows anonymized diagnostic patient samples to be used for purposes similar to those of the original sampling. The 29 samples from COVID-19 confirmed individuals, as well as the 100 negative controls and the 24 healthy volunteers were all from unique individuals. All samples were analyzed anonymously.

COVID-19 IgG/IgM rapid test

The test was run according to the manufacturers instructions (COVID-19 IgG/IgM Rapid Test Cassette (whole blood/serum/plasma), Product/Model: GCCOV-402a, Lot: 2003242, Zhejiang Orient Gene Biotech Co Ltd, Huzhou, Zhejiang, China) [11]. Briefly, 5 µl of serum or one drop of capillary blood were added to the test slide, followed by 80 µl of the buffer provided in the kit. The results were read after 10 min (max 15 min), by the naked eye. Only tests in which the control line changed its color were regarded as valid (3 out of 156 (1.9%) cassettes did not function). If a line was observed for IgM and/or IgG, the test was considered positive. The intensity of the color was not judged.

Results

IgM and IgG reactivities in negative control samples

None of the 80 negative sera from healthy blood donors tested IgM positive in the assay, while one tested IgG positive (1/80, 1.25%, 95% confidence level: 0.03–6.77%) (Tables 1 and 2). The single IgG-positive sample was re-analyzed and remained IgG positive in the second test. None of the 20 serum samples from the 6–12 months old babies tested positive for either IgM or IgG.

Table 1. Comparisons of IgM results for 29 PCR-positive COVID-19 cases and 124 healthy individuals.

	Cases	Healthy	Total
IgM positive	20	0	20
IgM negative	9	124	133
Total	29	124	153

Table 2. Comparisons of IgG results for 29 PCR-positive COVID-19 cases and 124 healthy individuals.

	Cases	Healthy	Total
IgG positive	27	1	28
IgG negative	2	123	125
Total	29	124	153

IgM and IgG reactivities in healthy volunteers

None of the 24 healthy volunteers, without any known history of SARS-CoV-2 infection/COVID-19, tested positive for IgM or IgG.

IgM and IgG reactivities in PCR-confirmed COVID-19 patients

Altogether 20 of 29 (69%) samples from PCR-confirmed COVID-19 patients tested IgM positive and 27 tested (93.1%) IgG positive (Tables 1 and 2). When the patients were grouped into two groups depending on the time between onset of disease and testing, seven out of ten patients in the first group (9–17 days) and 13/19 patients in the second group (18–29 days) tested IgM positive. Nine out of ten patients in the first group (9–17 days) and 18/19 patients in the second group (18–29 days) tested IgG positive (Figure 1). There was no statistical difference between the two groups for neither IgM or IgG seropositivity. All samples that were IgM positive were also IgG positive.

Assay sensitivity, specificity and accuracy

Based on the results described above and summarized in Tables 1 and 2, the assay showed a sensitivity of 69% (20/29) and 93.1% (27/29) for IgM and IgG, respectively. The assay showed an overall specificity of 100% (124/124) and 99.2% (123/124, 1 false positive) for IgM and IgG, respectively.

Using PCR-positive cases as true positives, the accuracy of the test was 94.1% (144/153) and 98.0% (150/153) for IgM and IgG, respectively. The positive and negative predictive values (the likelihood of being a case given a positive test result, and the likelihood of being healthy given a negative test result) for IgM were 100% (20/20) and 93.2% (124/133), respectively.

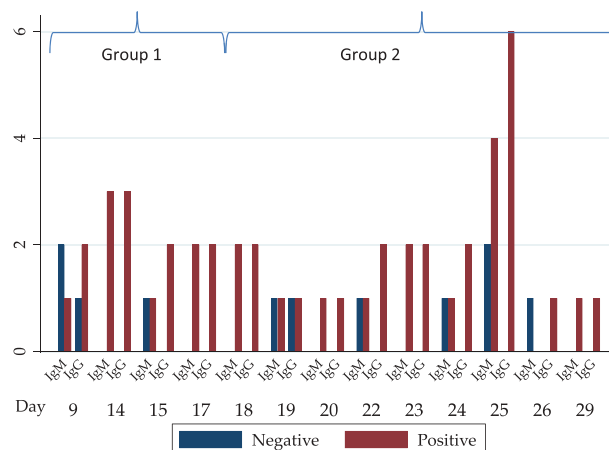


Figure 1. Number of PCR-positive cases positive or negative for IgM or IgG based on number of days after onset of COVID-19 symptoms.

For IgG, the corresponding values were 96.4% (27/28) and 98.4% (123/125).

Discussion

In this study we evaluated a commercial rapid test for detection of SARS-CoV-2-specific IgM and IgG. For the evaluation, samples from COVID-19 cases, obtained during disease or convalescence and previously confirmed by PCR, were used as 'true positives'. This means that in the PCR-positive cases for which antibodies may not yet had time to develop, or in potential cases with immune defects, it is possible that the negative IgM or IgG results were in fact true negatives. If this was the case for one or more of the included patients, the actual sensitivities should be higher, i.e. when evaluated only on samples known to contain detectable levels of SARS-CoV-2-specific IgM and/or IgG. For a more optimal evaluation of the assay sensitivity, a gold standard for SARS-CoV-2-specific antibodies would have been needed. This is, however, unfortunately not available to date.

According to the manufacturer, the specificity has been evaluated on 14 PCR-negative samples and was found to be 100% for both IgM and IgG, while the sensitivity evaluated on COVID-19 cases was calculated at 87.9% for IgM and 97.2% for IgG. The results by Li et al. [11] indicated an overall testing sensitivity of 88.7% and 90.6% specificity. Our results showed a lower sensitivity for IgM, a similar sensitivity for IgG, and specificities in between the results of the two evaluations.

A recent study on three Chinese COVID-19 cases found that seroconversions occurred between 7 and 12 days after the onset of symptoms [12]. However, larger studies on the detailed kinetics of the antibody responses (e.g. IgA, IgM, IgG, neutralizing antibodies) are now urgently needed for a better understanding of the dynamics of the immune response during COVID-19. The results of our study showed detectable IgM and IgG in some patients at day 9, while in other patients the seroconversion seems to occur later. The impact of early or late seroconversion on the case severity is not known, and must now be explored. Interestingly, there were no IgM positives that were not IgG positive. Generally, IgM is produced first, and later there is a switch towards IgG production [13], but studies on SARS-CoV suggest that IgM and IgG often develop around the same time [14,15]. Our results are in line with this (Figure 1), but more detailed studies on long-term sequential samples from patients are now needed. It may be worth looking into whether this is a problem with the test, or a constant finding within the immune response to SARS-CoV-2.

There were no false IgM positive samples, indicating a very high specificity of the test. One false positive IgG result was observed for one healthy adult blood donor. This sample was re-tested and the result was consistent,

indicating a cross-reaction to another coronavirus. Serological cross-reactions have earlier been observed between SARS-CoV and SARS-CoV-2 [16]. There are other human coronaviruses (NL63, OC43, 229E, and HKU1) that are globally endemic or epidemic [17], and it may be possible that this reaction represented a cross-reaction due to a previous infection with one of those. Human CoV NL63 has been shown to use the same receptor, angiotensin-converting enzyme 2 (ACE2), as SARS-CoV and SARS-CoV-2 [18], which may indicate potential cross-reactive epitopes. How common the CoVs are as causative agents for 'common colds' is not known in detail, but there has been estimates that up to 20% of cases could be caused by CoVs [19].

The specificity and sensitivity for IgG detection of the rapid test evaluated here is well in line with those of a recently reported enzyme-linked immunosorbent assay (ELISA), which had a specificity and a sensitivity of 97.5% [20].

While this study showed a satisfactory performance of the rapid test, it is limited by being compared only to clinical cases and PCR-positivity, and as a next step, it is necessary to compare this assay to other serological tests. In contrast to Li et al. [11], we found less indications for using this test for clinical diagnosis. Nevertheless, it might contribute to detecting potential asymptomatic infections as well as getting a notion of the magnitude of the spread in different geographical areas, which might be a key to taking the appropriate decisions and policies forward. The high negative predictive value indicates that the rapid test will be useful for detecting past infections and possible immunity, which may be crucial for restoring social functions after lockdown.

Acknowledgments

This work was supported by the Swedish Research Council (VR, grant numbers 2016-02596, 2017-05807 and 2018-02569). The rapid tests that have enabled this study were donated to us by the Swedish company Noviral AB (organization number: 559175-7942).

Disclosure statement

The authors declare no conflicts of interest.

Notes on contributors

Tove Hoffman is a PhD student at the Department of Medical Sciences, Uppsala University, Sweden. Her doctoral studies mainly concern the involvement of northward migratory birds in the dispersal of ticks and tick-borne pathogens.

Karolina Nissen is a medical resident and responsible for preventive measures for High Consequence Infectious

Diseases at the Infectious Disease Clinic, Uppsala University Hospital.

Janina Krambrich is a PhD student at the Department of Medical Biochemistry and Microbiology, Uppsala University, Sweden. Her doctoral studies concern molecular epidemiology and overwintering strategies of Japanese encephalitis virus.

Bengt Rönnberg, PhD, is working as a researcher at the Laboratory of Clinical Microbiology at Uppsala University Hospital and at the Zoonosis Science Center at Uppsala University. His main research focus is development of diagnostic assays and studies on zoonotic viruses.

Dario Akaberi is a PhD student at the Department of Medical Biochemistry and Microbiology, Uppsala University, Sweden. His doctoral studies mainly concern research on antivirals against flavi- and coronaviruses.

Mouna Esmaeilzadeh, MD, PhD, is a medical doctor and neuroscientist from the Karolinska Institutet, Sweden. She was one of the founders of SciLife Clinic (2009-2016), a pioneering medical clinic within preventive and personalized medicine, using technologies such as magnetic resonance imaging and genetic sequencing to map and prevent diseases. In addition to her MD, PhD, she has a master's degree in philosophy from the University of Oslo and Stockholm University.

Erik Salaneck is an associate professor in Infectious Diseases, Department of Medical Sciences, Uppsala University and responsible physician for High Consequence Infectious Diseases at the Infectious Diseases Clinic, Uppsala University Hospital.

Johanna Lindahl is a veterinary epidemiologist working at Zoonosis Science Centre at Uppsala University, with affiliations at International Livestock Research Institute (ILRI) in Nairobi Kenya, and Swedish University of Agricultural Sciences. Johanna graduated from Swedish University of Agricultural Sciences after doing her PhD working on Japanese encephalitis virus in Vietnam. Since her PhD she has been focusing her research on food safety, and vector-borne, zoonotic and emerging infectious diseases in developing countries, mainly East Africa and South and Southeast Asia.

Åke Lundkvist, PhD, is Professor in virology and one of the founders of the Zoonosis Science Center at Uppsala University. His work focuses on basic virological, pathogenetic, ecological and epidemiological studies on zoonotic viruses and their vectors and hosts.

ORCID

Tove Hoffman  <http://orcid.org/0000-0003-2239-9078>
 Janina Krambrich  <http://orcid.org/0000-0002-5061-7143>
 Bengt Rönnberg  <http://orcid.org/0000-0003-0229-093X>
 Erik Salaneck  <http://orcid.org/0000-0002-0000-8948>
 Johanna Lindahl  <http://orcid.org/0000-0002-1175-0398>
 Åke Lundkvist  <http://orcid.org/0000-0001-8608-6551>

References

- [1] Guarner J. Three emerging coronaviruses in two decades: the story of SARS, MERS, and Now COVID-19. *Am J Clin Pathol*. 2020 Mar 9;153(4):420–421.
- [2] De Wit E, Van Doremalen N, Falzarano D, et al. SARS and MERS: recent insights into emerging coronaviruses [Internet]. *Nature Reviews Microbiology* Nature Publishing Group, London, UK; Aug 1 2016;523–534. Available from: www.nature.com/nrmicro
- [3] Hilgenfeld R, Peiris M. From SARS to MERS: 10 years of research on highly pathogenic human coronaviruses. *Antiviral Res*. 2013;100:286–295. [Elsevier].
- [4] Wang D, Hu B, Hu C, et al. Clinical Characteristics of 138 hospitalized patients with 2019 novel coronavirus-infected pneumonia in Wuhan, China. *J Am Med Assoc*. 2020;323(11):1061. .
- [5] Huang C, Wang Y, Li X, et al. Clinical features of patients infected with 2019 novel coronavirus in Wuhan, China. *Lancet*. 2020 Feb 15;395(10223):497–506.
- [6] Wu Z, McGoogan JM. Characteristics of and important lessons from the coronavirus disease 2019 (COVID-19) outbreak in China: summary of a report of 72314 cases from the Chinese center for disease control and prevention. *J Am Med Assoc*. 2020;323(13):1239.
- [7] Li Q, Guan X, Wu P, et al. Early transmission dynamics in Wuhan, China, of novel coronavirus-infected pneumonia. *N Engl J Med*. 2020 Jan 29;382(13):1199–1207.
- [8] Wu X, Cai Y, Huang X, et al. Co-infection with SARS-CoV-2 and influenza A virus in patient with pneumonia, China. *Emerg Infect Dis*. 2020 Jun;26(6). doi:10.3201/eid2606.200299.
- [9] Worldometer. Worldometer, Coronavirus [Internet]. 2020 [cited 2020 Apr 4]. Available from: <https://www.worldometers.info/coronavirus/>
- [10] European Centre for Disease Prevention and Control (ECDC). Novel coronavirus disease 2019 (COVID-19) pandemic: increased transmission in the EU/EEA and the UK – sixth update [Internet]. Stockholm, Sweden; 2020. Available from: <https://www.ecdc.europa.eu/sites/default/files/documents/RRA-sixth-update-Outbreak-of-novel-coronavirus-disease-2019-COVID-19.pdf>
- [11] Li Z, Yi Y, Luo X, et al. Development and clinical application of a rapid IgM-IgG combined antibody test for SARS-CoV-2 infection diagnosis. *J Med Virol* [Internet]. 2020 Feb 27. cited 2020 Apr 4; jmv:25727. Available from: <https://onlinelibrary.wiley.com/doi/abs/10.1002/jmv.25727>
- [12] Zhang J, Liu J, Li N, et al. Serological detection of 2019-nCoV respond to the epidemic: A useful complement to nucleic acid testing. *medRxiv*. 2020;10:20030916.
- [13] Gunther N, Hoffmann GW. Qualitative dynamics of a network model of regulation of the immune system: A rationale for the IgM to IgG switch. *J Theor Biol*. 1982 Feb 21;94(4):815–855.
- [14] Li G, Chen X, Xu A. Profile of specific antibodies to the SARS-associated coronavirus. *N Engl J Med*. 2003;349:508–509.
- [15] Hsueh PR, Huang LM, Chen PJ, et al. Chronological evolution of IgM, IgA, IgG and neutralisation antibodies after infection with SARS-associated coronavirus. *Clin Microbiol Infect*. 2004 Dec 1;10(12):1062–1066.
- [16] Wan WY, Lim SH, Seng EH. Cross-reaction of sera from COVID-19 patients with SARS-CoV assays.

- medRxiv [Internet]. 2020 Mar 23. cited 2020 Apr 1;20034454. Available from: <https://www.medrxiv.org/content/10.1101/2020.03.17.20034454v1>
- [17] Paules CI, Marston HD, Fauci AS. Coronavirus infections-more than just the common cold. *J Am Med Assoc.* 2020;323:707–708.
- [18] Hofmann H, Pyrc K, Van Der Hoek L, et al. Human coronavirus NL63 employs the severe acute respiratory syndrome coronavirus receptor for cellular entry. *Proc Natl Acad Sci U S A.* 2005 May 31;102 (22):7988–7993.
- [19] Bende M, Barrow I, Heptonstall J, et al. Changes in human nasal mucosa during experimental coronavirus common colds. *Acta Otolaryngol.* 1989;107(3–4):262–269. .
- [20] Zhao R, Li M, Song H, et al. Serological diagnostic kit of SARS-CoV-2 antibodies using CHO-expressed full-length SARS-CoV-2 S1 proteins. *medRxiv.* 2020;27:20042184.